

# Drive X


Built to Perform Made to Last



Dental Implant Unit


Gullin Woodpecker Medical Instrument Co.,Ltd.  
Information Industrial Park, Gullin National High-Tech  
Zone, Gullin, Guangxi, 541004 P.R. China  
Sales Dept.: +86-0773-5873196/2350599  
After-sales Service Dept.: +86-0773-5827898  
E-mail: woodpecker4@gullinwoodpecker.com  
Website: http://www.gullinwoodpecker.com

DISTRIBUTED BY:



## Outstanding quality Smooth power

Advanced motor control technology blends low-speed precision with high-speed stability for confident and reliable performance

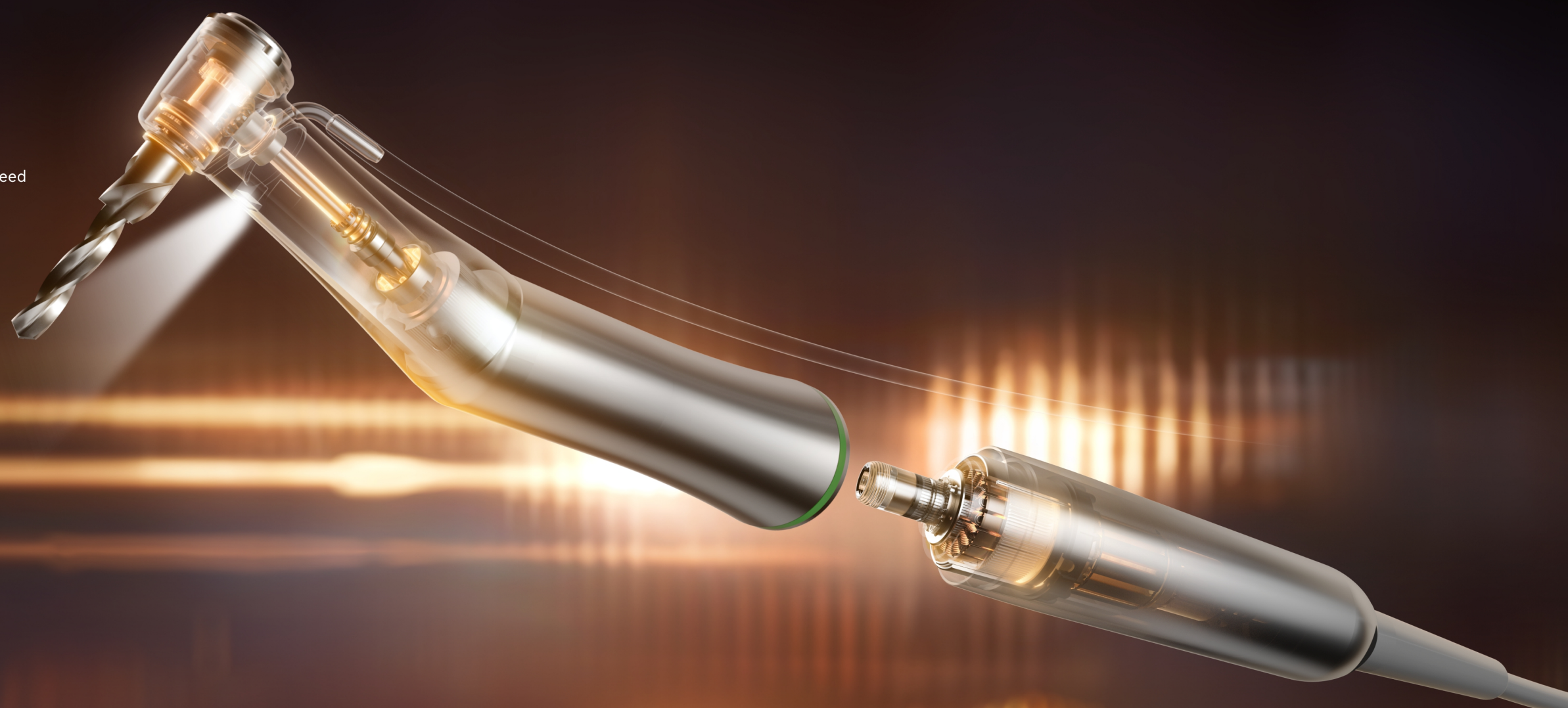

### LED contra-angle for clearer visibility

Built-in LED for clear illumination, enhancing visibility and precision during procedures




### Durable contra-angle with quality bearings

Fitted with 8 high-quality Chinese bearings for smooth operation and long-term durability under high use frequency




### Precision contra-angle for stable power and comfort

Engineered for accuracy and comfort — with steady power and minimal vibration

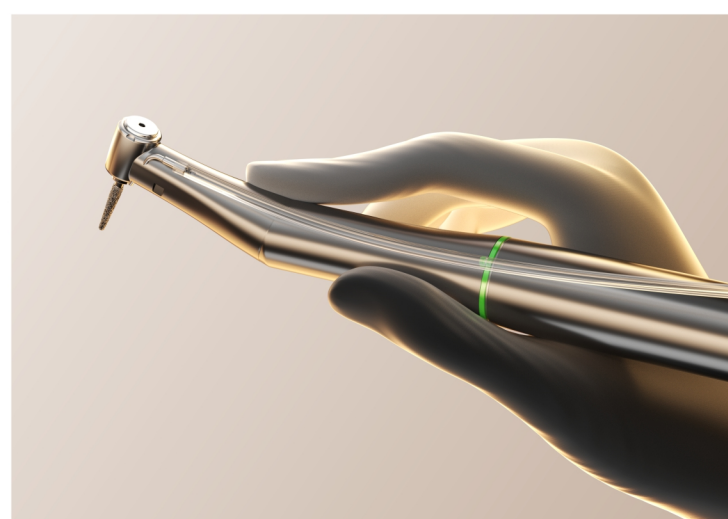

### High-performance domestic motor for confident precision

•6.2N·cm torque with up to 80N·cm output  
•Lightweight and autoclavable for smooth, safe implant procedures

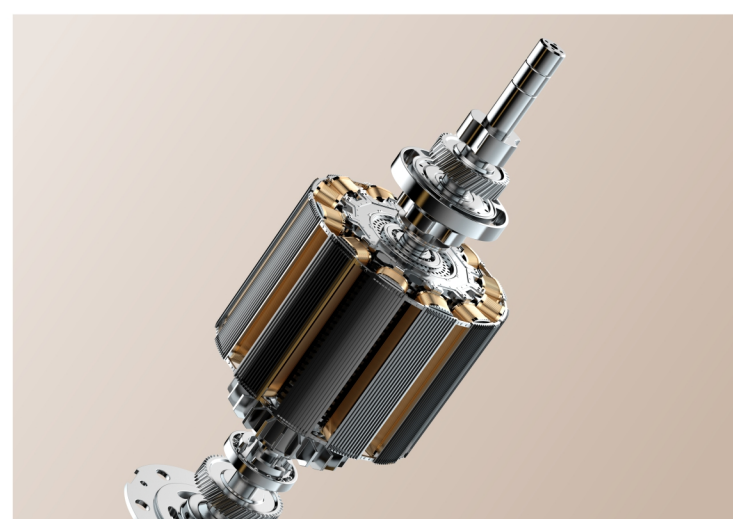

### Clear display Easy control

The large color screen shows every step clearly with preset settings ready to use

### Torque maintenance function

Stable torque control for precise and consistent implant procedures



### Brightness control

Switch between 50% and 100% LED brightness to suit different surgical lighting needs



### LED light delayed off

Set LED light delay to remain on briefly after operation, allowing time for post-procedure observation



### Volume control


Enable or disable sound prompts for a quieter clinical experience



Torque range/ Gear 20:1	5 N·cm – 80 N·cm
Max pump output	110mL/min
Main unit dimensions	250 mm * 237 mm * 130 mm
Motor dimensions	φ 23.4 mm * L 91 mm
Motor weight	150 g
Contra-angle head weight	88.3 g
Power supply	AC 220V, 50/60Hz

### Flexible foot pedal for hands-free operation

Hands-free control with flexible foot pedal for safer, more efficient procedures

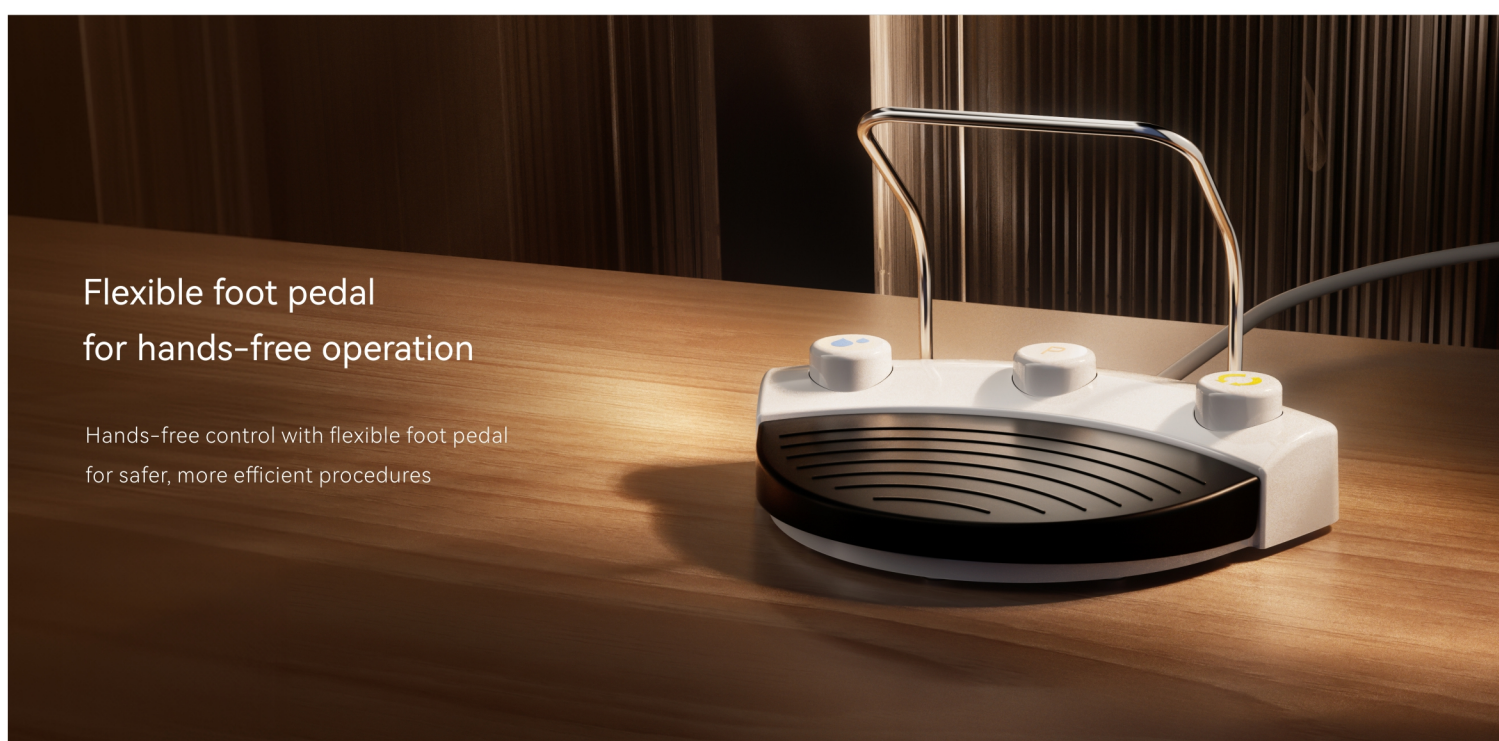

### Upgraded mute peristaltic pump

Featuring a four-roller design for more continuous and stable water delivery

